

Ultrasonic Homogenizer|Processor|Sonicator-1500W

Model No.: BEM-1500D

(LCD Screen Version)



Application

The instrument is widely used in the fields of biology, microbiology, physics, zoology, agronomy, pharmacy, petroleum, etc. Also suitable for sample preparation in the chemical, medical, archaeological, address, marine and other fields.

It also can be used in Nanotechnology research such as the dispersion of nanomaterials (nanotubes, graphene, silica, etc.); emulsification, homogenization and fragmentation of nano-materials such as carbon nanotubes, graphene and other materials, homogeneous soil, rock samples, non-ferrous metals, rare earths and other particles.Lysis, disruption, leaching, extraction of DNA, protein extraction from animal and plant cells, pruning DNA/RNA, chromatin immunoprecipitation technology for cells, viruses, spores and other cell structures.

Fax: +86- 25-51816989




Fax: +86- 25-51816989

Product Description

The Ultrasonic Homogenizer/Sonicator/Liquid Processor (LCD Screen) is a multifunctional and multipurpose instrument for ultrasonic processing by cavitation in liquid material. It can be used for a variety of plants and animals, bacteria, viruses, cells and the organization of broken, at the same time can be used to emulsify, separation, scattered, cracking, homogenization, extraction, defoaming, cleaning, nanometer material preparation, graphene dispersion and chemical reactions, etc.

Product Features

- 1. Large LCD screen Display, friendly human-machine interface.
- 2. **50 groups of experimental data** can be set and stored, safe and convenient.
- 3. The ultrasonic **power can be adjusted** continuously, and the best experimental conditions can be explored.
- 4. Efficient and durable titanium alloy probe, reliable and accurate.
- 5. The program has over-temperature/delay/fault protection settings and **alarm system,** which can protect the test samples to the maximum extent.
- 6. RS232 Interface is optional, which can support the data transmission between the PC and PLC to realize the remote control.
- 7. We provide the **Sound-proof Chamber for standard configuration** with perfect transportation.

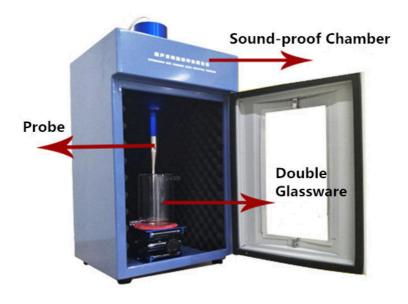



Fax: +86- 25-51816989

Technical Specification

Model	BEM-1500D		
Frequency (KHz)	19-21		
Power(Watt)	10-1500		
Processing Capacity (mI)	20-1000		
Input &Display	LCD Screen		
Login Password	Yes		
Standard Probe(mm)	Ф22		
Optional Probes (mm)	Φ10、15、20、25		
Sample Temperature Protection (°C)	0-100		
Security System	Self-diagnosis function, automatic error correction, overload, overtemperature protection, etc.		



Fax: +86- 25-51816989

Probe Options



Model	Inch	Frequency	Power Range	Crushing Capacity
Ф2	1/12"	20-25KHz	min-150W	0.2-5ml
Ф3	1/8"	20-25KHz	min-250W	3-10ml
Ф6	1/4"	20-25KHz	20-400W	10-100ml
Ф10	5/12"	20-25KHz	100-600W	30-300ml
Ф12	1/2"	20-25KHz	200-900W	50-500ml
Ф15	5/8"	20-25KHz	300-1000W	100-600ml
Ф20	4/5"	19.5-20.5KHz	400-1100W	100-1000ml
Ф22	5/6"	19.5-20.5KHz	400-1100W	200-1000ml
Ф25	1"	19.5-20.5KHz	800-1500W	500-1200ml